


mindray

BC-20s

Auto Hematology Analyzer

**Minimum Size
Maximum Capability**

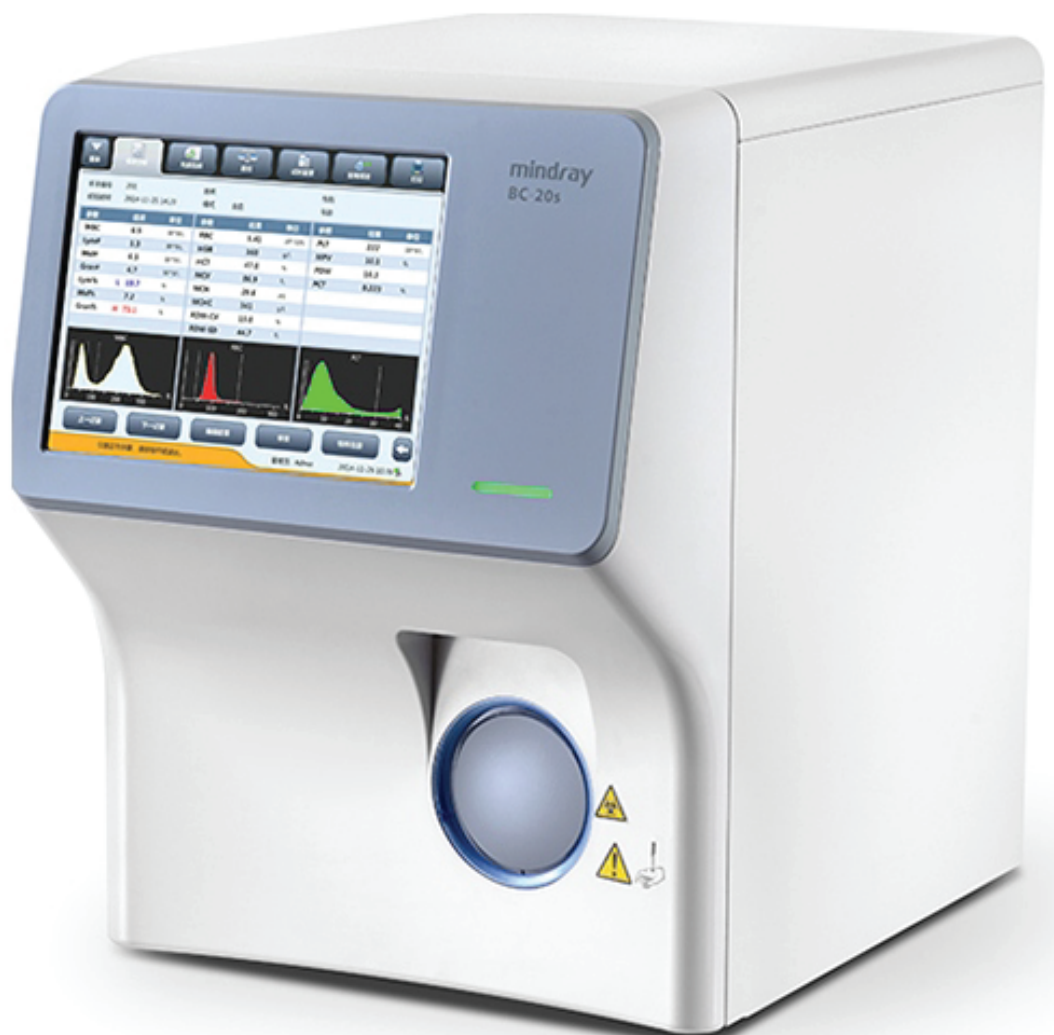

Technical Specifications

Principles

Impedance method for WBC,RBC and PLT counting
Cyanide free reagent for hemoglobin test

Performance

Parameter	Linearity Range	Precision (CV%)	Carryover
WBC($10^9/L$)	0-100	$\leq 2\%$ (7.0-15.0)	$\leq 0.5\%$
RBC($10^{12}/L$)	0-8.00	$\leq 1.5\%$ (3.5-6.5)	$\leq 0.5\%$
HGB(g/L)	0-280	$\leq 1.5\%$ (100-180)	$\leq 0.5\%$
PLT($10^9/L$)	0-1000	$\leq 4.0\%$ (150-500)	$\leq 1.0\%$

Parameters

19 parameters:WBC, Lymph#,Mid#,Gran#,Lymph%, Mid%, Gran%,RBC, HGB, HCT, MCV, MCH, MCHC, RDW-CV, RDW-SD,PLT, MPV, PDW, PCT
3 histograms for WBC, RBC and PLT
2 research parameters: NLR,PLR

Reagent

M-30D DILUENT
M-30CFL LYSE
PROBE CLEANSER

Sample Volume

Prediluted mode	20 μ L
Whole blood mode	9 μ L

Throughput

40 samples per hour

Display

8.4 inch TFT Touch Screen

Data Storage Capacity

Up to 200,000 results including numeric and graphical information

Communication

LAN Port supports HL7 protocol
Support bi-directional LIS

Interface

4 USBport (for external printer, software upgrade, barcode reader, external WIFI, keyboard and mouse), LAN port (1)

Printout

Thermal recorder, 50mm wide paper, various printouts formats
External Printer optional

Operating Environment

Temperature: 10°C~35°C
Humidity: 10%~90%
Air pressure : 70kPa ~106kPa

Power Requirement

100V-240V
 $\leq 300VA$
50Hz/60Hz

Dimension and Weight

Dimension: Depth(410 mm)x width(300 mm)x height(400 mm)
Weight: $\leq 20Kg$

UBM
UNITRON BIO-MEDICALS

UNITRON BIO-MEDICALS

AN ISO 9001:2015 CERTIFIED COMPANY

"Shree Durga", #464/1, 12th Main Road, 17th "A" Cross
JP Nagar 5th Phase, Bengaluru - 560078 India

Tel: 080-25187600
Customer Care: 9449009531
Email: info@unitronbiomedicals.com
web: www.unitronbiomedicals.com